

Model	Victory® 10 <sup>1</sup>
Product Number	3W: SC610 4W: SC710
Weight Capacity	400 lbs.
Maximum Speed <sup>2</sup>	Up to 5.3 mph
Ground Clearance <sup>4</sup>	3W: 2.2" at motor, 1.26" at anti-tip 3.4" at side deck 4W: 2" at motor, 1.5" at anti-tip 3.1" at side deck
Turning Radius <sup>4</sup>	3W: 42" 4W: 54"
Overall Length <sup>4</sup>	3W: 45.25" 4W: 47"
Overall Width <sup>4</sup>	3W" 22.25" 4W: 22.25"
Seat-to-Ground Height	23.5"-25"
Seat-to-Deck Height	17.5"-19"
Front Tires	3.5" x 10" solid
Rear Tires	3.5" x 10" solid
Suspension	N/A
Range Per Charge <sup>2,3</sup>	Up to 15.5 miles
Total Weight <sup>6</sup>	Without batteries: 3W 133 lbs.; 4W: 137 lbs. With batteries: 3W-181 lbs. (U1); 3W-197 lbs. (40AH); 4W-185 lbs. (U1); 4W-201 lbs. (40AH);
Heaviest Piece When Disassembled <sup>4</sup>	3W: 54 lbs. (front section) 4W: 63 lbs. (front section)
Standard Seating	Type: Compact, foldable Weight: 28 lbs. Material: Black Vinyl Dimensions: 18" width, 17" depth
Drive System	Rear-wheel drive sealed transaxle with 24-volt DC motor
Dual Braking System	Regenerative and electromechanical
Battery Requirements <sup>5,7</sup>	Size: (2) U1 standard, Weight: 24.5 lbs. each Size: (2) 40 AH (optional), Weight: 32.5 lbs. each
Battery Charger	3.5-amp, off-board (U1), 5-amp, off-board (40 AH)
Warranty	Frame: lifetime limited Drivetrain: 2-year limited Electronics: 2-year limited Battery: 6-months

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

**PRIDE**<sup>®</sup>  
MOBILITY  
LIVE YOUR BEST<sup>®</sup>

(US) 800-800-8586 • (Canada) 888-570-1113

pridemobility.com |    

HCPCS Code: K0801 (Group 1, Heavy Duty)

- 1 Pride® FDA Class II Medical Devices are designed to aid individuals with mobility impairments.
- 2 Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).
- 3 Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).
- 4 AGM or gel cell type required.
- 5 Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- 6 Battery weight may vary based on manufacturer.
- 7 Includes standard seat weight.

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

0003984-012722

**VICTORY**<sup>®</sup> 10

**ENHANCE  
YOUR ACTIVE  
LIFESTYLE  
WITH ADDED  
COMFORT AND  
PERFORMANCE**



FDA Class II Medical Device<sup>1</sup>

**LIVE YOUR BEST<sup>®</sup>**

pridemobility.com |    

# STYLE, COMFORT, PERFORMANCE

18" x 17"  
Pinchless hinge,  
black vinyl seat  
with 8-point  
swivel and  
sliders for  
increased  
comfort



Exclusive stylish, lightweight, non-scuffing, black, low-profile tires




Ambient voltmeter

Standard front basket

Easy-to-grip tiller adjustment knob

Long-lasting, maximum-intensity LED headlight offers optimal pathway illumination

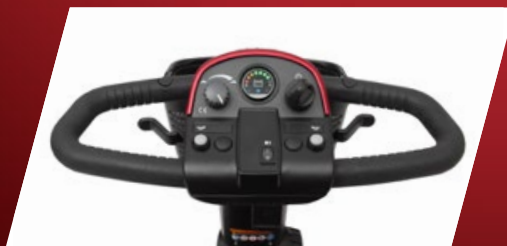

Tiller console



Feather-touch disassembly

The Victory® 10 series offers added style, comfort and performance for your active lifestyle. Standard features include LED pathway light, LED battery meter, feather touch disassembly, wraparound delta tiller, solid tires and a front basket. A 400 lb. weight capacity, a per charge range of up to 15.5 miles and a maximum speed up to 5.3 mph compliment the versatility of the Victory 10. Available in Viper Blue and Candy Apple Red.

- Exclusive one-hand feather touch disassembly
- Charger port located on the tiller lets you conveniently charge your scooter
- Delta tiller with ergonomic wraparound handles lets you operate the scooter with one hand and rest your wrist
- Easy-to-read ambient voltmeter
- Long-lasting, maximum-intensity LED headlight offers optimal pathway illumination

## Colors



Candy Apple Red



Viper Blue

## Optional Seating

Limited recline, black vinyl, high back with height adjustable headrest, flip-up height and width adjustable armrests

- 18" x 18" (sliders)
- 20" x 20" • 22" x 18"

Deluxe contour high back black vinyl with height adjustable headrest, flip-up height and width adjustable armrests (no sliders)

- 18" x 18" • 18" x 20"
- 20" x 20" • 22" x 20"

Pinchless black vinyl medium back hinge seat with 8-point swivel and sliders for increased comfort

- 20" x 18"

## Accessories

- Cane clip attachment
- Cane/crutch holder (single, double)
- Cell phone holder
- Cup holder
- Forearm crutch holder
- Oxygen tank holder
- Quad cane holder
- Rear basket
- Additional rearview mirror
- Saddle bag
- Safety flag
- Walker holder
- Weather cover
- XLR USB charger
- Lap belts: 50", 60" and 70"

## Military Patches

